



CALVING....



Rope Economy

single eye each 201 061 \$12.95
double eye each 201 060 \$19.95

Lowest cost calving rope option. 7mm polypropylene rope. Single version is 1.1m long with eye splice one end. Two of these should be used for assisting calving. The double version is 1.6m long with eye splices at both ends. Supplied each.

PRICE REDUCED



Rope Lightweight Pair

209 915 ~~\$32.95~~ \$29.95

Lightweight 4mm dia nylon ropes with eye splices both ends, sealed for durability. Length 1.5m o.a. Supplied as a pair (two ropes).



Rope Slimline

each 201 063 \$22.95

Soft synthetic rope, 8mm x 1.5m specially made for obstetrics. Small finely spliced eye loop both ends. Excellent for internal uterine use. Two required for correct use.



Rope Flat Braid

9mm 213 186 \$16.95
12mm 211 076 \$19.95

Soft, safe, flat-braid rope made especially for Shoof. Featuring 'stitched flat' areas before each eye loop. These create flat 'pads' around the ankles of the calf so as to prevent any tissue damage when pulling. Length is 1.8m for each, long enough to go around your back for easier pulling.




Strap Heavy Duty

each 201 064 \$19.95

All purpose calving strap of double-layer heavy-duty 25mm webbing with stainless steel D-rings double-stitched each end. Length 1.5m o.a. Best used as a pair.

CALVING CHAIN LENGTHS

After considerable testing we have chosen 80cm (short) and 150cm (long) as ideal calving/obstetric chain lengths. The short chain is intended for use in a pair, whilst the long chain can be used on its own, with each end connected to a leg. (We recommend the short chains, used as a pair, as the most appropriate for the job). When fitting chains to calve's legs it is best to use a double-loop to reduce possible injury.

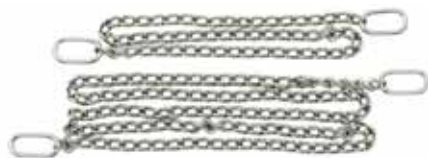

Chain Nickel-plated

short/ea 80cm 201 054 \$39.95
long/ea 150cm 201 053 \$59.95

Our nickel plated calving chains are of top quality USA manufacture, tested to over 500kg. Hardened, and nickel plating for long life. Oval end loop for easy internal use. Short chains should be used in pairs.



Chain Stainless Quality

short/ea 80cm 201 059 \$34.95
long/ea 150cm 201 058 \$59.95

These are our best quality chains. The 150cm chain features an intermediate loop-link for greater versatility especially when used with a pulling device. 160kg tested breaking strain.

WARNING ABOUT CALVING CHAINS

Calving chains should not be used with mechanical calf pullers, for the following reasons:

1. The chain links can do serious damage to the calf's leg.
2. The chains may break under the excessive force sometimes applied.

Use ropes with mechanical calving aids.

CHAIN CORROSION

All calving chains will corrode, including stainless steel chains, if left lying in antiseptic solutions for extended periods. To prevent corrosion, wash chain after use, and store dry.



Obstetric Gown Heavy Duty

204 880 \$199.00

No need to have to wash all your clothes after a calving job. Arms-free adjustable sleeve seals. Neck to boots with no openings. Waterproof and easily cleaned. Use for any obstetrical jobs. Long life rip-free laminated fabric in discreet colour. One size only, approximately medium to large.



Chain Handle

each 201 051 \$24.95

Strong stainless steel handle to use with our calving chains. Hooks in end loops or between any links. Top quality manufacture. Plastic grip. Two required for correct use.



Chain Handle (discontinued)

each 201 052 ~~\$59.95~~

Specially cast stainless steel handle with improved grip. Hook will attach in any chain link or loop. Top quality manufacture. Two required for correct use.



Strap Heavy Duty Comfort

each 201 065 \$24.95

The Comfort version of our calving strap has special additional hand-grip loops 50cm out from the ring ends. The strap is 2.5m long overall, allowing it to go around the back of the user. When assisting calving, especially when the cow is lying down, the user can use his or her legs, and the hand grips to assist purchase. Made from double-layered 25mm webbing with stainless steel D-rings.



Umbilical Cord Clamps

50-pack 210 697 \$15.00
100-pack 210 698 \$24.95

Research has confirmed that several diseases, including E.coli scours and Joint-ill, find their way into the newborn animal via the umbilical cord. Sprays and dips can also cause problems. Our umbilical clamps are the solution. Simply attach as soon as possible after birth. Remove after 24 to 48 hours. Place clamp approximately 5mm from top of umbilical cord. Cut cord below clamp to prevent suckling. Clamps can be disinfected and re-used if desired, or can be left to fall off naturally. Ideal for piglets as clamps are large enough not to fall through creep floors. Clamps 65mm long o.a.

Copyright © 2011. All rights reserved.